




# LASER THERAPY (SCANNING and PORTABLE)

Model Digi Laser - Pro  
 - Gives Guaranteed Result

Universal Power Supply: 110 - 275 VAC 50Hz / 60Hz  
 Fuse Rating :1A ±10%  
 Max. Power Consumption : 30W  
 Safety Standard IEC60601 - 1  
 Safetyclass Class 1 type BF (IEC60601 - 1)  
 Electrical and functional safety test certified

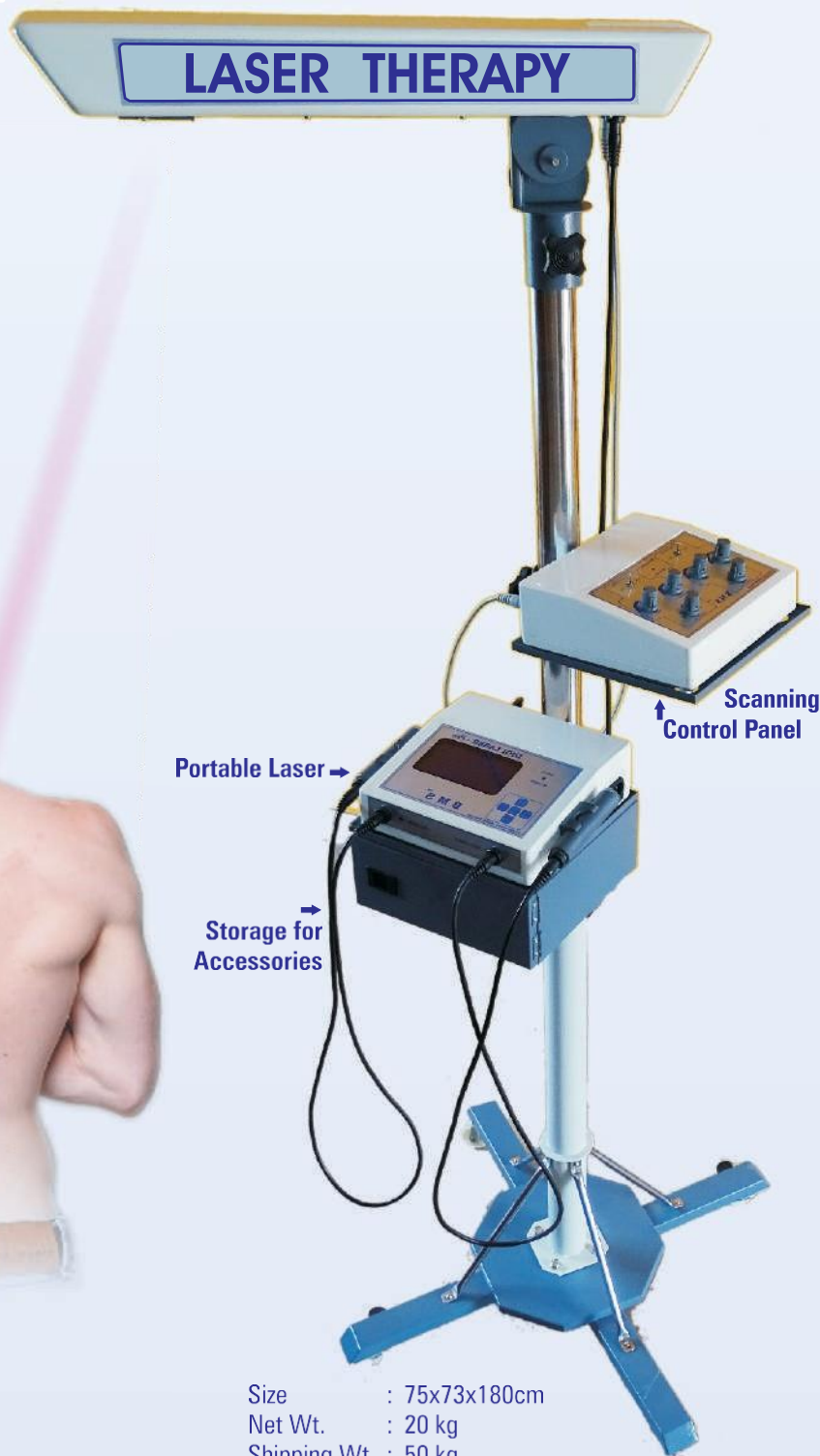

Size : 75x73x180cm  
 Net Wt. : 20 kg  
 Shipping Wt. : 50 kg

## SPECIFICATION

This can be used as Scanning Laser Therapy as well as Portable Laser Therapy.  
 Arm can be tilt up 0-60° and Down 0-60° by Manual. (Remote adjustments is optional) Arm can be rotate 360°  
 Scanning Laser Control Unit and Portable Laser Therapy Unit Height can be adjusted, storage box for accessories.  
 computerized Control panel can set the various pattern of Laser beam shape and size upto X' axis 2 feet and Y' axis 2 feet Independently.  
 Laser Beam Pattern like circle, ellipse, Horizontal, Vertical, Spot and combination.  
 Scanning area minimum 1 cm<sup>2</sup> to 360 cm<sup>2</sup>, Scan the set area automatically and delivers the even dosage of energy to the treatment area, it provides faster healing.